
 SILVERGEAR

SMART SCALE PRO

BODY COMPOSITION & HEART RATE

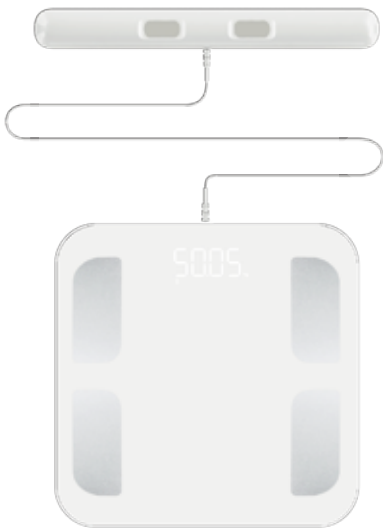**User Manual
Handleiding**

EN	2
NL	10
FR	18
DE	26

Thank you for purchasing our smart scale with handle bar. This manual will provide instructions for effective and safe use of the product. Please read it carefully before use and store the manual in a safe place for future reference. If you ever pass the product on to others, this manual should be included.

Table of contents

- 1.Package content
- 2.Technical data
- 3.Description of parts
- 4.Safety instructons
- 5.Operation
- 6.Maintenance
- 7.Troubleshooting
- 8.Disposal

1. Package content


- Scale
- Handle bar
- Manual

2. Technical data

Product size:	scale: 300 x 300 x 27.4 mm handle bar: 280 x 34 x 32 mm display: 83 x 35 mm
Weight unit:	Kg/lb/st:lb
Weight limit:	180 kg/400 lb/28 st:lb
Weight division:	0.1 kg/0.2 lb
Power supply:	4x 1.5V AAA LR03 batteries (not incl.)
Operating temperature:	10-40 °C
Operation humidity:	40-80 % RH
Support devices:	iOS: 8.0 or higher Android: 6.0 or higher
Maximum Radiated Power:	-1.4 dBm
Operating Frequency:	2402MHz-2480MHz

3. Description of parts

1. Electrode sensors
2. Measuring handle bar
3. Connecting wire
4. LED display
5. Sensor feet
6. Battery cover

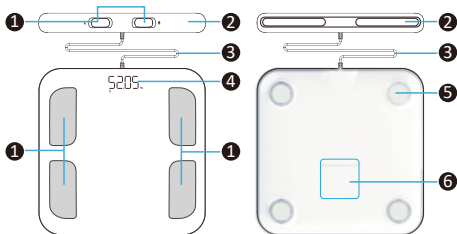4. Safety instructions

- The product is safe to use by children from the age of 10 and by persons with reduced physical or mental capabilities, or with lack of experience or knowledge, but only under supervision of an adult and with taking into account any hazards.
- Maintenance should be done by or under the supervision of an adult.
- Before use, inspect the device and its parts. Do not use in case of damage.
- Do not use to diagnose or treat any medical conditions.
- Do not use with medical implants such as pacemakers.
- Keep the cord out of harm's way.
- Make sure the scale and your feet are dry before stepping on the scale, to avoid slipping.
- Do not stand on the edge of the scale or jump on it.
- Do not overload the scale (max 180 kg/400 lb).
- Place the product on a hard level surface.
- Do not drop the scale or drop objects on it, as this may damage the sensors.
- Do not immerse in water.
- Not recommended for pregnant women.
- Not recommended for infants, toddlers, and children under 10 years of age.
- Do not attempt to repair, disassemble or alter the product yourself.

Note

- Any measurement obtained using this device, is for reference only and should not be considered as a medical opinion.
- Consult your GP or specialist before making changes to your diet, exercise plan or physical activities.
- To ensure the most accurate results, take your measurements at the same time of the day and under similar circumstances.

5. Before use

1. Search for 'Fitdays' in the Apple Store or in Google Play. Or scan this QR code:

2. Download and install the app. The app icon will appear on your phone or tablet after the installation is completed.

3. Turn on Bluetooth and your GPS on your smart phone or tablet.


4. Register your own ID by email or login with a social media account.

5. Add personal data and confirm OK.

Note: Entering the wrong gender, age, height, or corporeality level can result in inaccurate measurements.

- Activate the scale by tapping on it with your foot. Wait until the display shows "0.0 kg".
- Pair the scale with your device via Bluetooth. Leave the app on the main screen, it will automatically connect to the scale.
Or click "Account"→"Device"→+→"Bluetooth" to find the scale's Bluetooth ID: My_scale.
- The app will show your data within a few seconds.

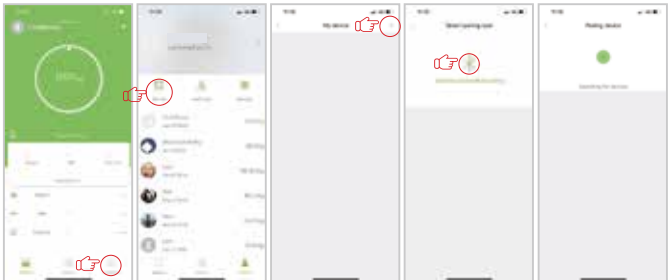6. Operation


Before using the smart scale, make sure you have entered all your personal data through the Fitdays app.

1. Place the scale on a hard level surface.
2. Tap with your foot to activate it. Wait until the display shows '0.0' (kg/lb).
3. Step on the scale barefoot for thorough measurements and full analysis.
4. Hold the handle bar in front of you with your arms straight, at a 45 degree angle. Make sure your arms do not touch your body, as this will negatively affect the measurements. Also, your hands and feet should make full contact with all sensors.
5. After the weight has been locked, all body composition data can be read in the Fitdays App by clicking 'more data'.

STEP 1

STEP 4


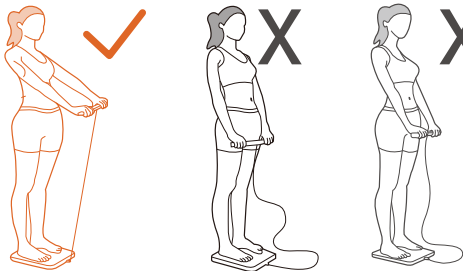
X

X

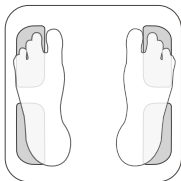X

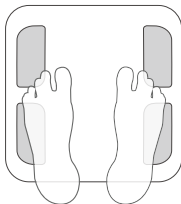X

Note: When stepping on the scale with socks or shoes on, only weight and BMI will be measured and displayed on the main page in the app.

7. Maintenance

- If necessary, clean the scale and handle bar with a slightly damp towel and wipe it dry with a soft cloth.
- You can use cleaning alcohol to clean the sensors. Do not use it on other parts of the unit.
- Store the scale in a dry, cool place away from direct sunlight. Be careful not to kink the connection wire.

8. Troubleshooting

1. Bluetooth connection failed.

- Make sure the scale, Bluetooth and your GPS are on.
- Check the system version of your smart phone, it must be iOS 8.0 or higher/ Android 6.0 higher.
- When failed to connect via Bluetooth, please click "Account" → "Device" → + → "Bluetooth" to find the scale's Bluetooth ID: My_scale.

2.The scale is not measuring body fat data.

- You have to be barefoot to measure body fat. Make sure your feet and the scale are both dry.
- Stand on the scale for about 10 seconds, until the number on the display stops flashing.
- Bluetooth is not turned on.
- Bluetooth is paired with a difference scale.

3.The scale is producing inconsistent measurements.

- Check whether the scale is placed on a hard, flat surface.
- Make sure nothing is stuck to the bottom.
- Remove the batteries to restart the scale.

4.The app is displaying an incorrect body fat rate.

- The body parameter settings are incorrect, check to make sure you entered the correct gender, height, and age.

5.How to re-calibrate the scale.

- Remove one of the batteries, place it back and wait for the scale to restart automatically.
- Tap with your foot on the scale to turn it on. Let it automatically drop to 0.0kg to calibrate.

6.The scale does not turn on.

- Remove the batteries and place them back. Make sure to insert the batteries as directed.
- Check if the batteries are depleted.

For more information, please check FAQ from the Fitdays app.

LO: The batteries in the scales are too weak. Insert new batteries.

EER: The scale is overloaded (max 180 kg). Place less weight onto the scale.

C: An error had occurred during measurement. Step off the scale for 5 seconds and step back on to measure again.

9. Disposal

Dispose of the article and the packaging materials in accordance with current local regulations.

The crossed out wheeled bin symbol means that this product shall not be disposed of with normal household waste. Electronic and Electrical Equipment not included in the selective sorting process are potentially dangerous for the environment and human health due to the presence of hazardous substances. Please dispose of responsibly at an approved waste or recycling facility.

The CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. The CE marking is the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

Recycling helps to reduce the consumption of resources and reduce the environmental impact. If you wish to discard the packaging, comply to the environmental legislation in your country.

Simplified Declaration of Conformity

- Hereby, Karsten International declares that Item number: 3052, Smart Scale with handle bar, is in compliance with the RED Directive 2014/53/EU. The full text of the EU DoC is available at the following internet address:

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