130mm

HARDCOVER FRONT AND BACK

 SILVERGEAR

MASSAGE GUN PREMIUM

EN - USER MANUAL NL - HANDLEIDING FR - GUIDE DE L'UTILISATEUR DE - BEDIENUNGSANLEITUNG

www.silvergear.eu

Overschiestraat 63
1062 XD Amsterdam
The Netherlands
info@silvergear.eu

Not For Resale
Made in China

Art. no.: 2723
Batch no.: 21939

130mm

Thank you for purchasing our Massage gun premium. This manual will provide instructions for effective and safe use of the product. Please read it carefully before use and store the manual in a safe place for future reference. If you ever pass the product on to others, this manual should be included.

Table of contents

1. Package content
2. Technical data
3. Description of parts
4. Safety instructions
5. Operation
6. Maintenance
7. Disposal

1. Package content




- Massage gun
- Adapter
- 6 Massage heads
- Storage case
- Manual

2. Technical data


Colour: black
 Material: PS, ABS
 Size: 17x5x24cm
 Weight: 1.1kg
 Input: 100-240V AC, 50/60Hz 5A
 Output: 26V DC 1A
 Battery capacity: 2400 mAh
 Charging time: 2.4 hours
 Working time: about 5 hours (on level 1)
 Speed range: 1600-3200 rpm (30 levels)
 Low noise level: 55dB
 Smart timer: automatic shut-off after 10 minutes

3. Description of parts

1. Motor
2. Battery
3. LED light
4. Charging port
5. Power switch
6. Touchscreen
7. Speed level
8. Battery level
9. Power button
10. Speed +
11. Speed -

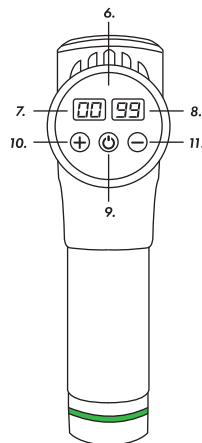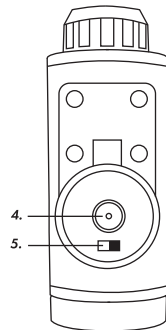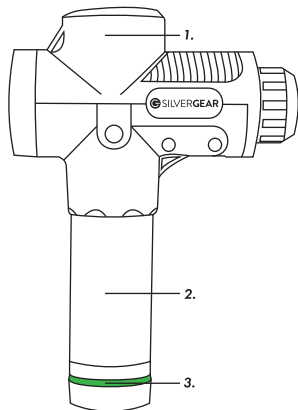Massage heads

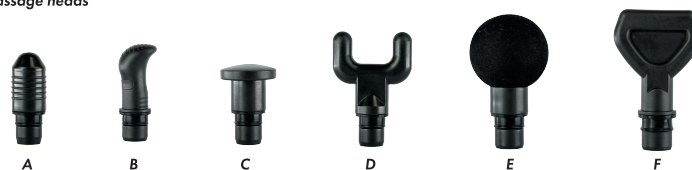- A.** Spherical head: for large muscle groups, such as arms, thighs, etc.
B. Thumb head: for soft tissue massage.
C. Flat head: to relax the muscles.
D. U-Shape head: for neck, spine and Achilles tendon.
E. Cylindrical head: for deep tissue massage of joint.
F. Shovel head: for the back (cervical, lumbar massage).

4. Safety instructions

- For intended use only, as described in this manual.
- The product is safe to use by children from the age of 8 and by persons with reduced physical or mental capabilities, or with lack of experience or knowledge, but only under supervision of an adult and with taking into account any hazards.
- Cleaning and maintenance should be done by or under the supervision of an adult.
- Before use, inspect the device and its parts. Do not use in case of damage.
- Before use, consult your specialist in case of a medical condition or treatment.
- Not suitable for diabetics or people with poor blood circulation.
- Do not place near heat sources such as radiators, cookers, etc.
- Keep away from sources of ignition.
- Do not place, store or use the product near water.
- Do not immerse in water.
- Do not let moisture come into contact with any electrical parts.
- Never touch the device or its plug with wet hands.
- Stop use immediately and contact your specialist, if you experience any discomfort or pain.
- Do not use on sensitive, irritated or damaged skin.
- Keep hands away from moving parts.
- Remove jewellery and tie hair back to prevent getting caught.
- Do not insert any foreign objects into the ventilation opening of the device.
- Always unplug before installing, removing or adjusting any of the accessories and for long-term storage.
- Attachments not recommended by the manufacturer may cause hazard.
- Only use on soft tissues of the body.
- Misuse may result in injury, and damage to the device.
- Do not attempt to repair, disassemble or alter the product yourself.

5. Operation

Charging

Fully charge the battery before first use. This will take about 2.4 hours.

1. Plug the provided adapter into the massage gun and connect it to the mains.
2. The LED light on the handle will slowly flash green to indicate it is charging.
3. Once the battery is fully charged, the LED stops flashing and lights up green.

Note:

- Charge the massage gun at any time and at any battery level.
- It is not recommended to completely deplete the battery.
- When the battery level is less than 20%, the colour of the LED on the handle changes to red.

How to operate

Before use, make sure the massage gun is unplugged.

1. Attach a massage head by simply pushing it in.
2. Flip the power switch on the bottom to turn on the massage gun. All LEDs turn on.
3. Press  on the touchscreen to activate it. The default speed level is level 1 (1600 rpm).
4. Repeatedly press or hold  to increase or decrease the speed level.
5. To turn off the massage gun, press  and flip the power switch again.

Note:

- Only use directly on clean, dry skin. Do not use through clothes.
- Limit use to 60 seconds per area.
- Some bruising may occur regardless of control setting or applied pressure. Frequently check the treated area during use.

6. Maintenance

Make sure the power switch is in the OFF position. Clean the massage gun with a slightly damp towel and wipe it dry with a soft cloth. Store it in the supplied storage case and keep it a dry, cool place, away from direct sunlight.

7. Disposal

The crossed out wheeled bin symbol on this product is to remind you that electrical products should not be disposed of with household waste. Please dispose of it responsibly at an approved waste or recycling facility. Check with your local authority for recycling advice.

Dispose of all packaging materials in an environmentally friendly manner.

The CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. The CE marking is the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

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