# ***User Manual***

**Electric Multifunction Massager - U**

# User Manual

## Electric Multifunction Massager - U

Please read this user manual carefully to fully understand the usage of the Electric Multifunction Massager - U. Keep the manual for future reference since it contains important information.

### IMPORTANT SAFETY INSTRUCTIONS


This is an electronic product. To reduce the risk of electrocution, please follow these instructions below:

1. Always unplug or switch the Massager to the OFF position after using or before cleaning.
2. Never leave the Massager unattended when plugged in. Unplug from the outlet when not in use.
3. Do not operate the Massager under a blanket or other cover. Excessive heat can build up and cause burns or fire.
4. Do not allow children to use or operate.
5. Never operate the Massager if any part of the massager or plug has been damaged.
6. Do not operate the Massager after it has been dropped or damaged.
7. Do not pull the power cord to move the Massager or use the power cord as a grip.
8. Keep the power cord away from heated surfaces.
9. Never operate the Massager around liquid or wet areas.
10. For adult use only. Do not use massager in close proximity to loose clothing or jewelry. Keep long hair away from massager while in use.

**IMPORTANT:** Any individual who may be pregnant or recently gave birth, has a heart condition or a pacemaker, suffers from diabetes or other sensory impairment, phlebitis and/or thrombosis, joint dysfunction, inflammation or injuries, weak bones, osteoporosis, a high fever, an abnormal or curved spine, is at an increased risk of developing blood clots, has pins/screws/artificial joints or other medical devices implanted in his/her body, is being treated for a medical condition or physical injury or is under a doctor care should consult with a medical professional before using a massaging device designed for home use.


### USING YOUR MASSAGER

1. Place the Massager over your shoulders, ensuring placement is comfortable and the massaging components are over your shoulders.
2. Insert the Power Adapter into the DC wire connection outlet on the Massager.
3. Plug the power adapter into your household electrical plug.

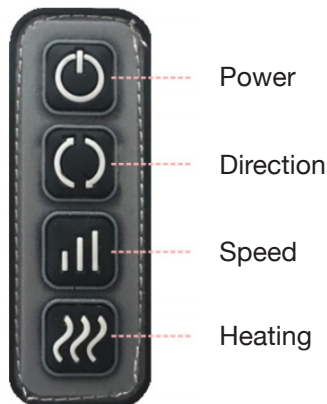
4. Turn the Massager ON by pressing the POWER button [1]. This will begin the shoulder massage program.
5. To make the massage nodes rotate clockwise, or counter-clockwise press the DIRECTION [2] button. Pressing the DIRECTION button will always change the direction of the massage nodes.
6. Change the speed of the massager and press [3] button. There are 3 gears in the speed adjustment, and the speed changes always change the speed.
7. To activate the HEAT option, press the HEAT button [4]. To deactivate the HEAT option, press the HEAT button again.

**NOTE:** The Massager includes a built-in 15 minute timer that activates every time the unit is turned on. This is the recommended time for massage use. If the unit turns off due to the automatic timer, and you wish to continue to use the Massager, simply turn the unit ON again.

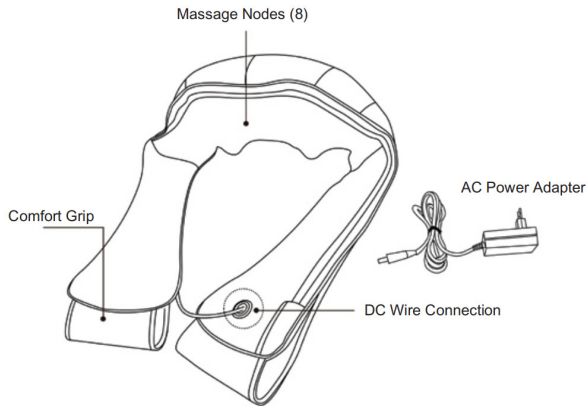**WHERE TO USE YOUR MASSAGER**

Although typically the Massager is to be used on your shoulders, it can also be used in many other spots on your body that are sore and require relief.
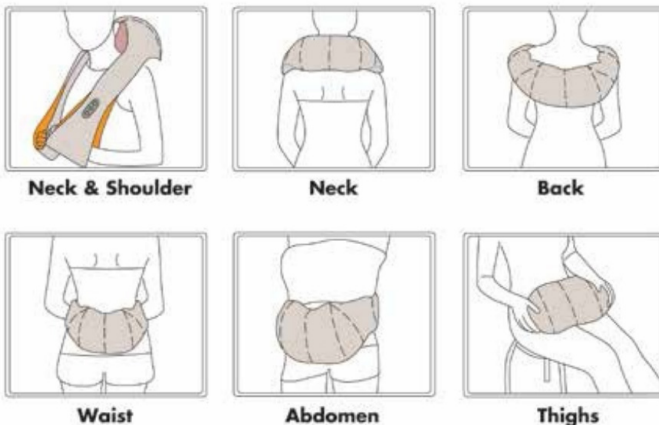
## **MAINTENANCE**

1. Before cleaning the product, always ensure the power plug has been removed from the unit and your wall outlet.
2. Use a lightly moistened cloth to wipe the unit off. Do not use abrasive cleaners on the Massager as this can cause unit damage or fading.
3. Never dunk or submerge the Massager in any liquids. Do not wash in the washing machine or other washing device.

## **TECHNISCHE SPECIFICATIES**

Package Size:	40 x 18 x 18cm
Massager Color:	Black
Gross Weight:	1.95 kg
Net Weight:	1.295 kg
Materials:	PU leather and wear-resistant breathable mesh
Rated Voltage:	DC 12V
Rated Power:	24W
Rated Frequency:	50Hz

Imported/Distributed by The Investment Company Europe B.V. exclusively developed for sales in the Netherlands, Belgium, and Luxembourg.

## **PATENT PENDING**

The artwork and design of this package are protected by European copyright law and may not be reproduced, distributed, displayed, published or used for any purpose without prior written permission. It is not permitted to alter, or reproduce any trademarks or copyright notice from this package.

NL / EN - [www.alleenonline.nl](http://www.alleenonline.nl)  
BE / FR - [www.alleenonline.be](http://www.alleenonline.be)

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