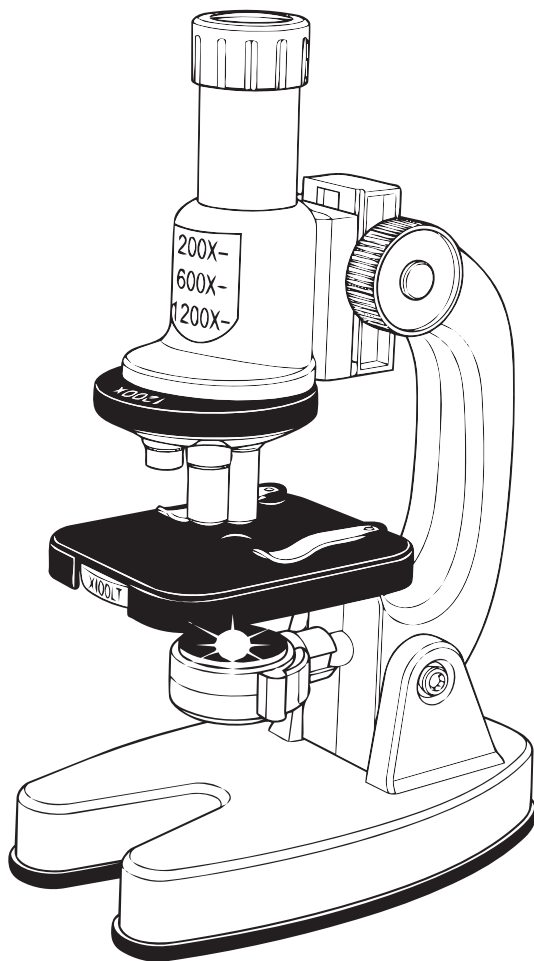
Microscope set

NL • Microscoopset DE • Mikroskop-Set
FR • Ensemble de microscopes EN • Microscope Set

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Important instructions. Retain for future reference.

Warnings


- Read the instructions before use, follow them and keep them for reference.
- Remove all packaging materials before giving the product to a child.
- WARNING:** Not suitable for children under 3 years. Small parts. Choking hazard.
- Non-rechargeable batteries are not to be recharged.
- Rechargeable batteries are to be removed from the product before being charged.
- Rechargeable batteries are only to be charged under adult supervision.
- Different types of batteries or new and used batteries are not to be mixed.
- Batteries are to be inserted with the correct polarity (+ and -).
- Exhausted batteries are to be removed from the product.
- The supply terminals are not to be short-circuited.
- Do not expose cells or batteries to heat or fire. Avoid storage in direct sunlight.
- Batteries can leak or explode if damaged, incorrectly installed, or exposed to heat.
- If the device is not used for a long time, remove the batteries.
- Dispose of the device, battery and all related components properly in accordance with local regulations.
- Check the product for damage before use. Do not use if damaged.
- Do not attempt any repairs or modifications; this may cause malfunction or injury.

Contents

- Microscope
- 20 Accessories
- Manual

Description of parts

1. Eyepiece
2. Lens barrel
3. Objective lens rotation
4. Objective lens
5. Stage
6. 4-colour filter
7. Light unit
8. Knob
9. Arm
10. Specimen clamp
11. ON/OFF switch
12. Base
13. Phone Holder

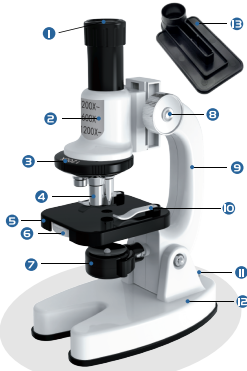
Description of accessories

1. Phone holder (1 PCS) – Holds a phone for taking pictures through the microscope.
2. Magnifying glass (1 PCS) – Enlarges objects for easier viewing.
3. Blank glass slides (10 PCS) – Hold samples for microscope viewing.
4. Duck feather glass slide (1 PCS) – Prepared slide for observing duck feather fibres.
5. Cotton glass slide (1 PCS) – Prepared slide for observing cotton fibres.
6. Coverslips (12 PCS) – Cover and protect samples on glass slides.
7. Concavity slide (1 PCS) – Hold small amounts of liquid for observation.
8. Slides storage box (1 PCS) – Stores and organises glass slides.
9. Coverslip storage box (1 PCS) – Stores and organises coverslips.
10. Dropper (1 PCS) – Adds small drops of liquid to samples.
11. Test tubes (2 PCS) – Hold small liquids or samples.
12. Petri dish (1 PCS) – Holds and observes small samples.
13. Specimen container (1 PCS) – Stores and transports small specimens.
14. Black cap small bottle (3 PCS) – Stores small samples or liquids.
15. Tweezers (1 PCS) – Picks up small objects safely and precisely.
16. Specimen knife (1 PCS) – Cuts small samples for observation.
17. Specimen shovel (1 PCS) – Collects small samples.
18. Specimen cone (1 PCS) – Moves or points to tiny samples.
19. Dust removal cloth (1 PCS) – Cleans the microscope and slides.
20. Anti-slip pad (1 PCS) – Keeps the microscope steady on the table.

Instructions

1. Insert 2 x 1.5V AA batteries into the base of the microscope and turn on the LED light using the ON/OFF switch.

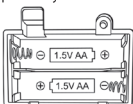2. If the LED light is too bright, adjust the position of the light unit. You can also use the mirror on the back of the light unit (no batteries required) by tilting and rotating it to reflect external light through the hole in the centre of the stage.
3. Place the specimen on the stage, secure it with the specimen clamps, and ensure the non-slip pads are attached.
4. Turn the focus knob until the objective lenses are at their highest position, then start observing with the lowest magnification lens and gradually move to higher magnifications up to 1200x. The higher the magnification of the objective lens, the greater the detail observed.
5. To focus the specimen, turn the knob until the objective lens is close to the specimen, then look through the eyepiece and slowly turn the knob in the opposite direction until the image becomes clear and fully focused.
6. Rotate the 4-colour filter located under the stage to switch between colours and experiment with different viewing effects.

Installation of phone holder

1. Remove the protective film from one side of the transparent stickers and attach them to the phone holder.
2. Place the phone holder on the eyepiece and remove the remaining protective film.
3. Turn on your phone camera, align the camera with the centre of the eyepiece lens, and place the phone securely on the holder.
4. Adjust the focus knob until the specimen image is clear on your phone screen, then take photos.

Maintenance

Clean lenses with lens paper only. Do not touch lenses with fingers. Keep the microscope dry and store it in a dust-free place. Handle with care.

Specifications

Zoom	200x, 600x, 1200x
Weight	0.38 kg
Age	6+
Batteries	2x 1.5V "AA" LR6 (not included)

Disposal

Dispose of the item and its packaging materials according to current local regulations.

This product is labelled with this crossed-out wheeled bin symbol in accordance with European Directive 2012/19/EC to indicate that it must not be disposed of with your other household waste. Due to the presence of hazardous substances, mixtures or components, electrical and electronic devices that are not subject to selective sorting are potentially dangerous to the environment and human health. Please check your local city office or waste disposal service for the return and recycling of this product.

The CE marking is a certification mark that indicates conformity with health, safety and environmental protection standards for products sold within the European Union. It represents the manufacturer's declaration that the product meets the requirements of the applicable EC directives.

Manufacturer

Overschiestraat 63
1062 XD, Amsterdam
The Netherlands
info@silvergear.eu
www.silvergear.eu

A digital version of this manual is available at:
<https://www.qrsten.com/>.

Overschiestraat 63
1062 XD Amsterdam
The Netherlands
info@karsten.nl

Item: 5559/3226307
Batch: 555926F24
Made in China